



26" HD Medical Display

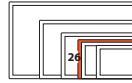


FSN Medical Technologies has solutions for managing your operating room video signals. Our products are engineered and built for compatibility with other highly specialized surgical and diagnostic equipment used in surgical suites, operating rooms, emergency rooms, and procedural facilities. FS-P2602D, our enhanced 26" medical display, is uniquely equipped to handle tasks in the demanding OR environment.

- Rapid signal detection, robust mode tables
- Artifact-free HD images
- Fanless - sterile field compatible
- Built to last, rugged aluminum housing
- Calibrated to clinical color
- Image pan, zoom, freeze, picture-in-picture

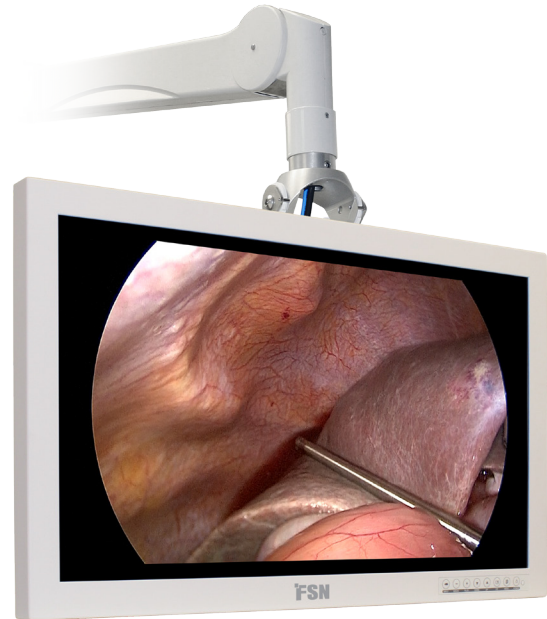


FSN has developed a process, seeking to improve picture quality, and applied it to our FS-P2602D displays. These enhanced LCD panels with LED backlight show an improvement in contrast ratio and brightness, plus improved humidity resistance and breakage resistance.

Accurate Color Reproduction

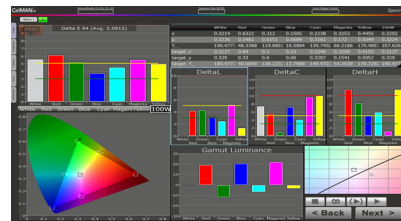


FSN has implemented methods to fine-tune the properties of our displays. Distinct color space settings have been calibrated to well-known surgical color preferences.

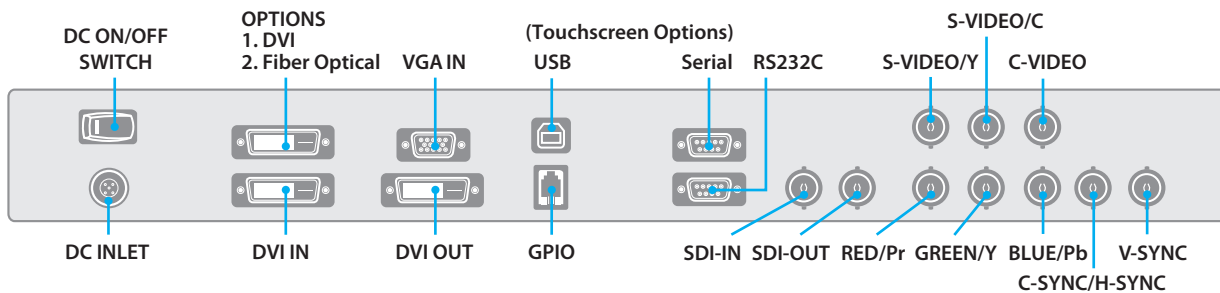

FSN displays have the connections you need. All major inputs are accepted via the back side connection panel.

Specifications

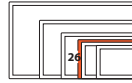


General Features

Item	Description
Panel	26 inch TFT LCD (IPS)
Resolution	1920 x 1080 pixel (HD 1080)
Aspect Ratio	16 : 9
Active Area	576.0 X 324.0 mm
Input Signal	1 x DVI-D 1 1 x DVI-D 2 (optical fiber optional) 1 x VGA (D-sub) 1 x SD/HD/3G-SDI (BNC) 1 x C-Video (BNC) 2 x S-Video (Y/C) (BNC) 1 x Component (RGBS, YPbPr) (5 x BNC)
TOUCH (optional)	(Note: Depth and weight will increase)
GPIO port	Swap, PIP / PBP1 / PBP2 Select, Record indicator
Output Signal	1 x DVI-D, 1 x SD/HD/3G-SDI (BNC)
Power Supply	AC/DC Adaptor (AC 100~240V, DC 24V 6.25A)
Control Key	Power, Menu, PIP, Up, Down, Plus, Minus, Input
Unit Dimension	638(W) x 389(H) x 74.7(D) mm 25.12(W) x 15.32(H) x 2.94(D) inch
Average Brightness	500 (5 points avg.) cd/m2
Response Time (G-to-G)	8 ms (avg.)
Viewing Angle (CR>10)	R/L 178, U/D 178
Weight	7.6 kg, 16.8 lbs
Contrast Ratio	1500 : 1 typical
Pixel Pitch	0.3 x 0.3

Signal Timing

Input Signal	Type	Description
VGA	Type	Analog RGB
	Connector	DSUB-15
	Level	0.7Vp-p ±5%
	Polarity	Positive
	Impedance	75 ohm ±5%
	Horizontal frequency	30 ~ 93KHz
	Vertical frequency	50 ~ 85Hz
	Applicable maximum pixel frequency	170MHz
	Type	Separate H/V sync, Composite H/V sync, Sync On Green Automatic synchronization without external switch
	Terminal resistance	more than 2KΩ
DVI	Type	Digital RGB
	Signal Format	TMDS single link
Optical DVI	Connector	DVI-D
	Type	Digital RGB
CVBS	Signal Format	Fiber Optic Transmission
	Type	C-Video
	Color system	NTSC, PAL
SVHS Y/C	Signal Format	Composite Video Signal
	Type	S-Video
	Color system	NTSC, PAL
SDI	Signal Format	Y/C Separated Video Signal
	SMPTE-425	1080p (60 / 50) 1080i (60 / 59.94 / 50)
	SMPTE-274M	1080p (30 / 29.97 / 25 / 24 / 24sF / 23.98 / 23.98sF)
	SMPTE-296M	720p (60 / 59.94 / 50)
	SMPTE-125M	480i (59.94)
ITU-R BT.656	576i (50)	

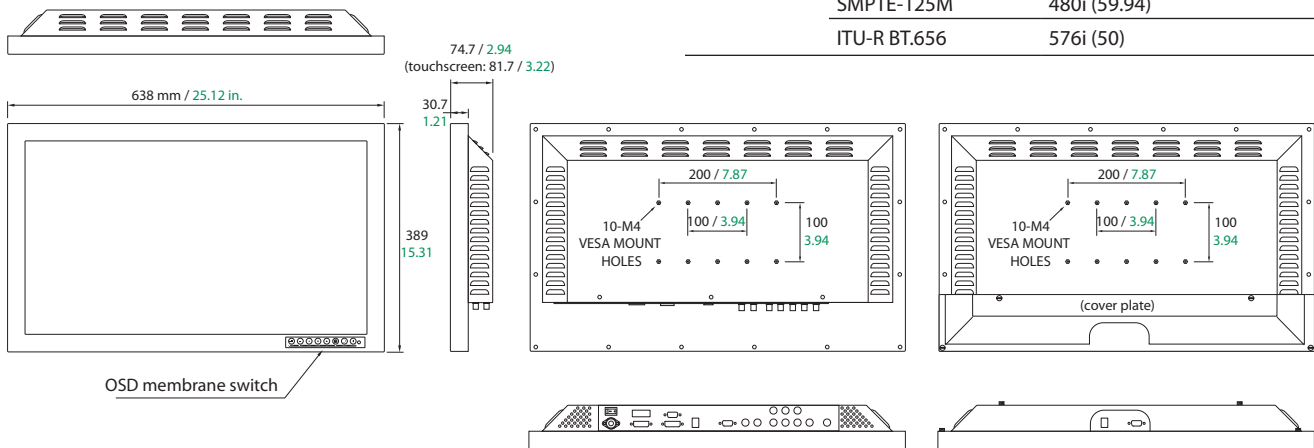

Specifications are subject to change with or without notice. Doc.# FSN1936 Rev. 08/14

www.fsnmed.com

Foreseeson Custom Displays, Inc.
2210 E. Winston Road
Anaheim, CA 92806 USA
Tel: 714-300-0540
Fax: 714-300-0546

Fort Lee, New Jersey
Tel: 201-849-4495
Atlanta, Georgia
Tel: 770-871-1895
Orlando, Florida
Tel: 407-667-3586

Foreseeson Korea
404B, Pangyo-innovaleys B
253 Pangyo-ro
Bundang-gu Seongnam-si
Gyeonggi-do, Korea 463-400
Tel: +82-31-8018-0780
Fax: +82-31-8018-0786

Foreseeson GmbH
Industriestrasse 38a
63150 Heusenstamm, Germany
Tel: +49 6104 64398 0
Fax: +49 6104 64398 11

Foreseeson UK Ltd.
Unit 2 Kingsmill Business Park
Chapel Mill Road
Kingston upon Thames, Surrey KT1 3GZ
Tel: +44 (0) 208 546 1047
Fax: +44 (0) 208 546 3931